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6 in 1 multi beauty device, aimed to address a wide range of skin concerns from A to Z for modern skincare

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Multifunction applicators provide skincare solutions, based on different skin types, tailored to various customer's needs, ensuring not only temporary improvements but fundamental skin enhancement and easy maintenance.



Product details	Description
핸드피스 Handpiece	SONOCLEAN, SONO(3way), HIFU(3way), RF, CRYO(COOLING&HEATING), ION ROD
본체 사이즈 Size	430(W) x 170(D) x 270(H)mm
본체 무게 Weight	10.0kg
정격전압 Rated voltage	AC 220V 150W / 48V3A 12V10A
소비전력 Power consumption	90W
LCD패널 LCD panel	800 x 480px (7inch)

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
The ultrasonic energy of 1, 3, 10, and 17MHz freq., vibrating at a rate of 1M times per sec., penetrates deep into the skin, enhancing elasticity without damaging. It stimulates collagen production, helping to create firmer and healthier skin. Additionally, it aids in the absorption of skincare products (such as ampoules), delivering moisture deep into the skin to make it radiance. Furthermore, it soothes the skin and promotes regeneration, aiding in the recovery of damaged skin, while activating lymphatic circulation to facilitate the efficient elimination of toxins, making it highly effective in improving problematic skin conditions.

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
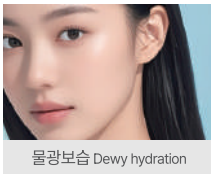
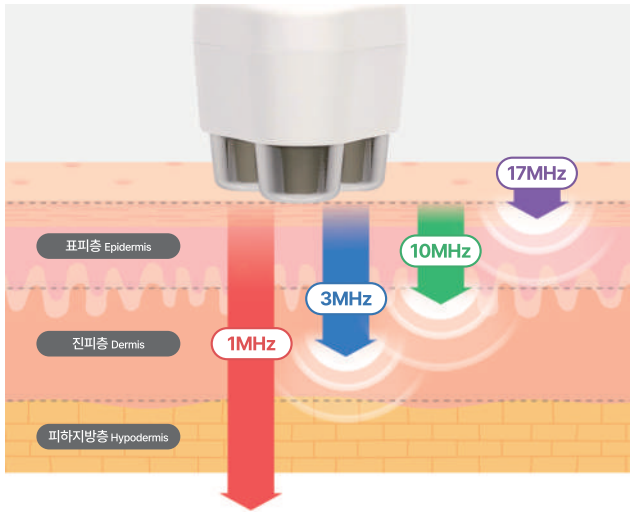
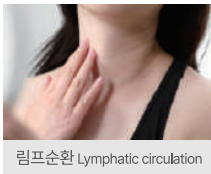
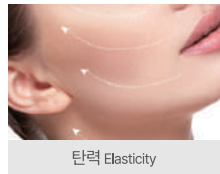
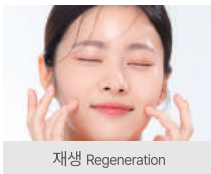
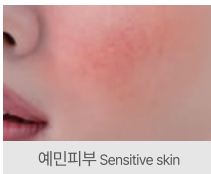
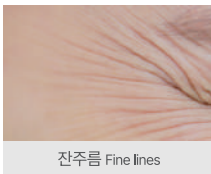
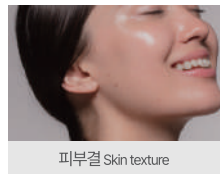
High-intensity focused ultrasound (HIFU) energy penetrates deep into the skin, contracting the SMAS layer (approximately 4.5 to 5mm deep) and providing an instant lifting effect. By effectively tightening sagging skin, it refines loosened facial contours and restores face definition, creating a polished facial outline. Over time, the skin keeps gradually tighten as the treatment continues to work beneath the skin surface, effectively reducing double chin and lifting drooping cheeks for a more youthful look. It provides non-invasive, firm lifting effect that fills in lost volume of the skin, smoothen wrinkles and defines contour, enhancing overall facial structure.

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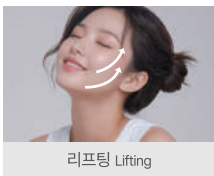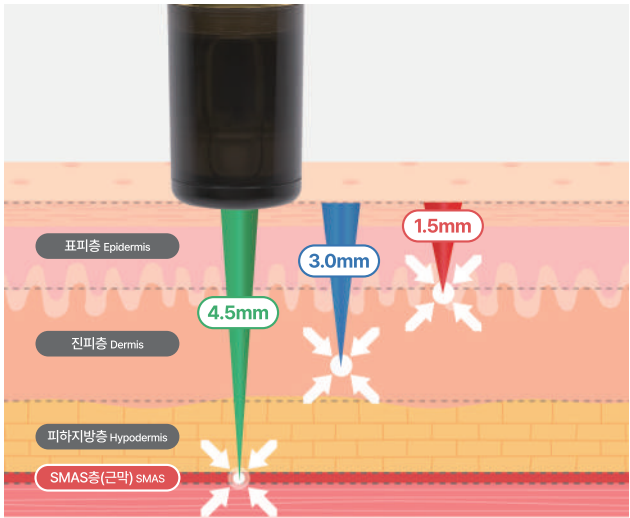
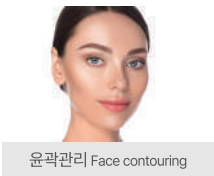
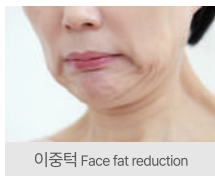
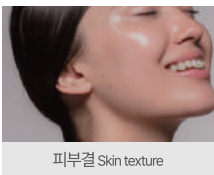
DOT thermal energy is delivered deep into the skin through 44 micro pins, providing intense thermo stimulation to activate collagen and elastin production, enhancing skin volume and elasticity. This process improves wrinkles, increases skin density, and aids in repairing damaged tissue, effectively reducing scars. Additionally, it tightens pores for smoother skin texture while promoting blood circulation, supporting cell regeneration, and revitalizing the skin. With consistent treatments, RF helps maintain firm, healthy, and radiant skin.

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



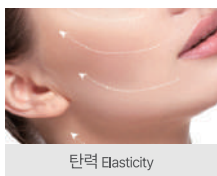
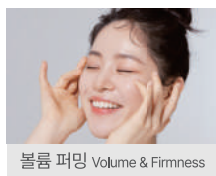
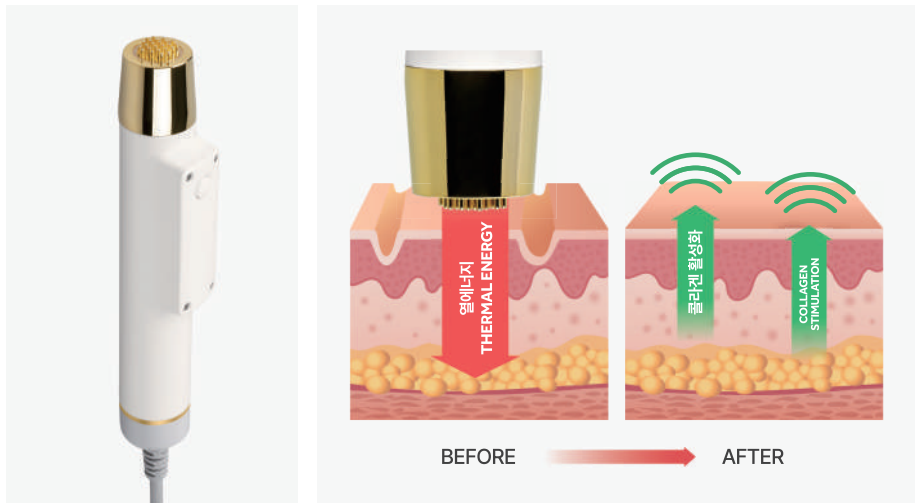
Using ultrasonic waves vibrating over 30,000 times per second at 34KHz, it effectively removes dead skin cells from the surface while deep clearing pores from sebum and blackheads to prevent acne formation. The micro-vibrations stimulate the skin, promoting blood circulation and boosting elastin and collagen production to improve skin elasticity and overall health. Additionally, it tones the skin, enhancing active ingredients penetration to the deeper layers of the skin. With regular use, it helps achieve a smooth and even skin texture and tone.

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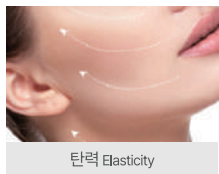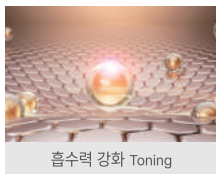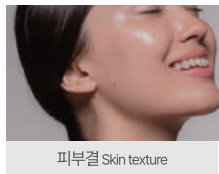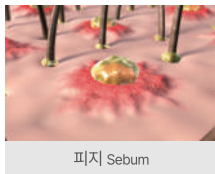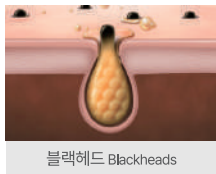
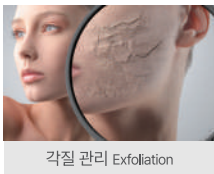
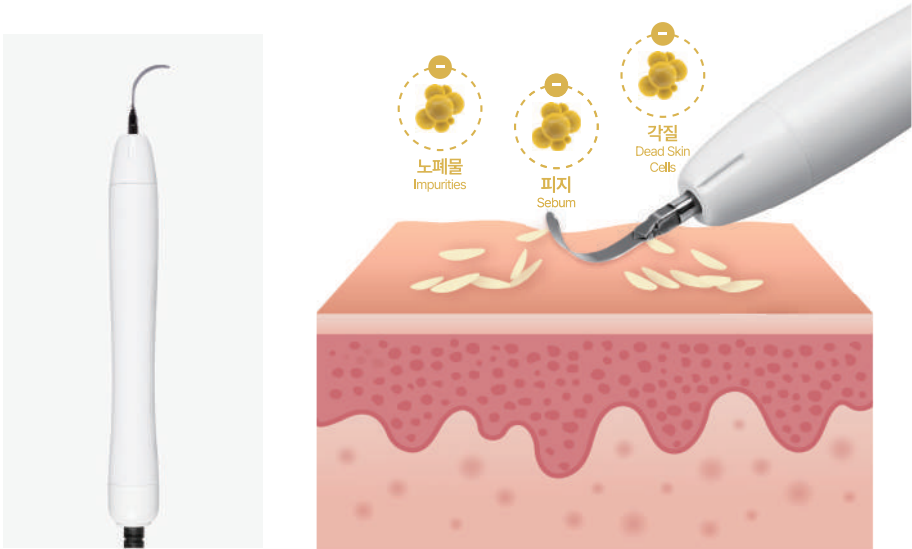
The cooling & heating therapy provides optimal skincare by adjusting the temperature according to the skin's condition, utilizing the Peltier effect based on thermoelectric principles. It is effective in soothing skin that has become sensitive due to laser treatments, microneedling, exfoliation, or exposure to UV rays and heat. Additionally, it can be used to open up the pores before acne extraction, making the process smoother and more efficient.

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


Electroporation is a skincare technology that utilizes micro-electric pulses to create temporary microscopic channels in the cell membrane, allowing active ingredients to penetrate deep into the skin. This technology effectively delivers high-molecular-weight substances such as hyaluronic acid and collagen, significantly enhancing absorption compared to conventional application methods. As a result, it maximizes the benefits of brightening skin tone, hydration, elasticity improvement, and post-acne scar reduction.

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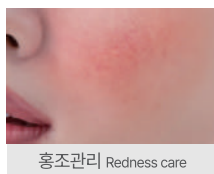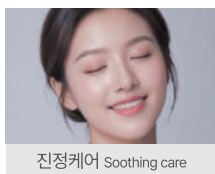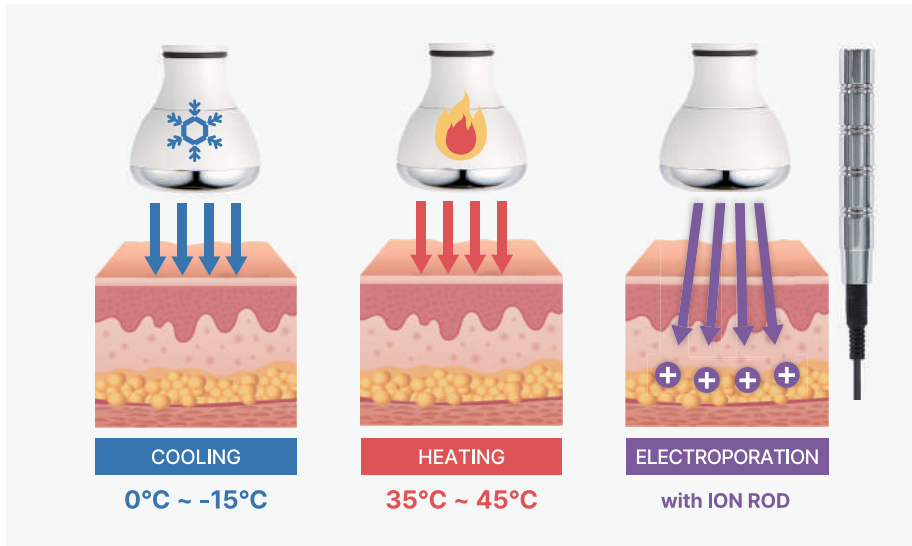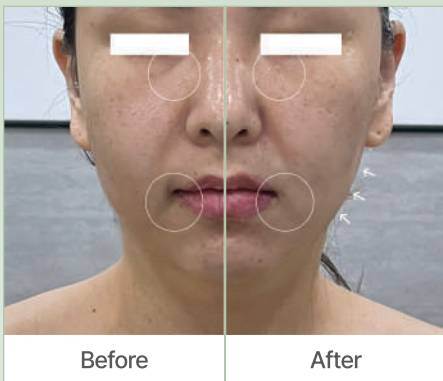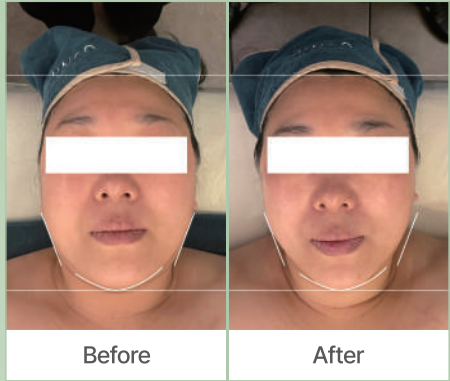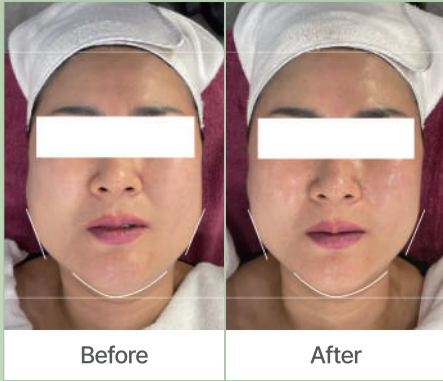
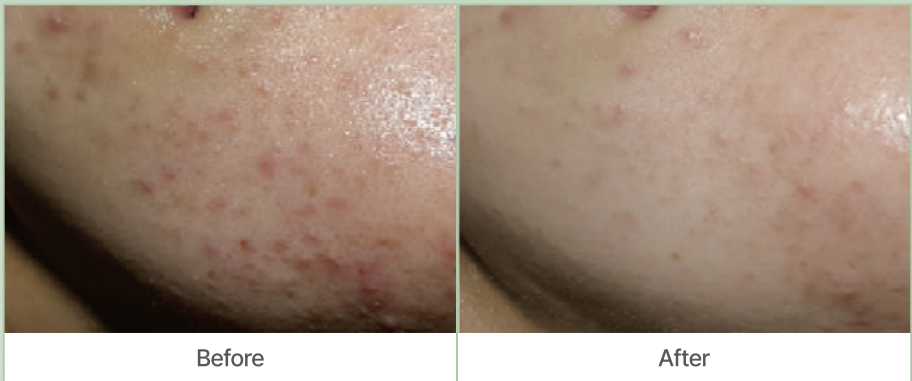
Precision Manufacturing. Proven Quality.

At Blue Solutions Medical Korea, every device is engineered through a meticulous manufacturing process that prioritizes precision, safety, and consistency. From high-quality component sourcing to advanced assembly and rigorous quality inspection, each stage is carefully controlled to meet professional standards.

Manufactured in Korea, our systems reflect our commitment to innovation, reliability, and excellence - ensuring optimal performance for clinics and aesthetic professionals worldwide.

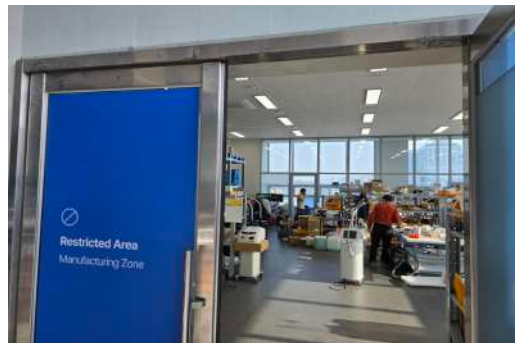BLUE SOLUTIONS MEDICAL KOREA

Your Vision, Our Solution.

Leading the Innovation in Beauty, Blue Solutions Medical Korea.

With cutting-edge technology and premium solutions,
we transform the future of aesthetic care.

Designed & Manufactured in Korea

📍 10, Jeongseojin 4-ro, Seo-gu, Incheon, Republic of Korea

☎ +82 10-8760-9812

🌐 www.bluesolutions.co.kr

✉ success@bluesolutions.co.kr